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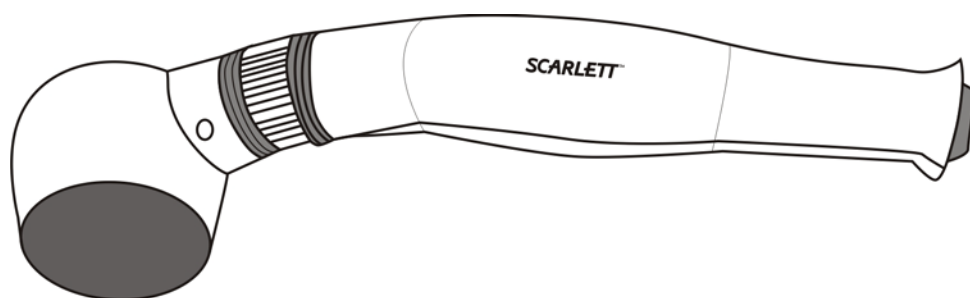
GB INSTRUCTION MANUAL

IMPORTANT SAFEGUARDS

- Please read this instruction manual carefully before use and keep in a safe place for future reference.
- Before the first connecting of the appliance check that voltage indicated on the rating label corresponds to the mains voltage in your home.
- For home use only. Do not use for industrial purposes. Do not use the appliance for any other purposes than described in this instruction manual.
- Do not use outdoors.
- Always unplug the appliance from the power supply before cleaning and when not in use.
- To prevent risk of electric shock and fire, do not immerse the appliance in water or any other liquids. If it has happened DO NOT TOUCH the appliance, unplug it immediately and check in a service center.
- Do not take this product to a bathroom or use near water.
- This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliances by a person responsible for their safety.
- Children should be supervised to ensure that they do not play with the appliance.
- Do not leave the appliance switched on when not in use.
- Do not use other attachments than those supplied.
- Do not operate after malfunction or cord damage.
- Do not attempt to repair, adjust or replace parts in the appliance. Check and repair the malfunctioning appliance in the nearest service center only.
- Keep the cord away from sharp edges and hot surfaces.
- Do not place or store the appliance where it can fall or be pulled into a tub or sink.
- Do not pull, twist, or wrap the power cord around the appliance.
- Never drop or insert any object into openings on the appliance.
- Do not put operating appliance on any surface, it may cause injury or damage.
- The massager is designed for intermittent use: it can be applied without interruption not longer than **20** minutes, then the motor must be allowed to cool.
- Do not use the massager if you are pregnant, or suffer from thrombosis, open wounds, broken skin and varicose veins. Consult your doctor first.
- Immediately interrupt procedure if you feel discomfort, or if any pain or irritation results. First, obligatory consult your physician.
- The appliance has a hot surface. Insensitive to heat persons should use the appliance carefully.
- If the product has been exposed to temperatures below 0°C for some time it should be kept at room temperature for at least 2 hours before turning it on.
- The manufacturer reserves the right to introduce minor changes into the product design without prior notice, unless such changes influence significantly the product safety, performance, and functions.

GENERAL INSTRUCTIONS

- Use of this massager has a relaxing, soothing and stimulating effect on your muscles. If used regularly for several minutes a day, it can reduce stress, fatigue and muscular pain, and it has anti-cellulite function.
- When using skin or face remedies, this appliance helps to spread the product more evenly and thus enhance its effective skin penetration.
- This appliance is equipped with an infrared lamp. The infrared radiation penetrates into the deep layers of skin causing the pores and vessels to expand and simultaneously enhancing the metabolism process and microcirculation. Using the infrared massager brings relief from muscular pain and regenerates the skin.

ACCESSORIES

ATTACHMENT FOR FACIAL MASSAGE

- This attachment is used for relaxation and toning facial muscles, increases the elasticity and resilience of the skin. Massage is performed in a circular motion. Avoid areas around the eyes and lips.

ATTACHMENT FOR ACUPRESSURE MASSAGE

- Acupressure massage helps relieve pain in muscles and tissues. Attachment for acupressure massage affects biologically active points of the body. It is recommended after heavy physical exercise.

ATTACHMENT FOR MICROPOINT MASSAGE

- Provides a more intense massage. Everyday use of this attachment will allow you to maintain high muscle tone and always be in good shape. Massage is performed in a circular motion. It is recommended to gently apply pressure to the skin to achieve better effect.

ATTACHMENT FOR ACUPUNCTURE MASSAGE

- This attachment is designed to gently massage large areas of the body (e.g., back, hips, etc.).

NOTE:

- In order to enhance the effectiveness of the massager move it slowly and evenly, in a smooth circular or back-and-forth motion.
- When massaging the body always start with those areas that are farthest from the heart, gradually bringing the device closer so that the heart is the center and the endpoint of the procedure.
- Never massage one part of the body continuously for more than three minutes.

ADJUSTMENT

- The angle of inclination of the massage head can be changed with the help of the rotating ring. Changing the angle you should watch the direction of rotation as shown in the figure and not overtwist the ring in order not to damage the device.

INSTRUCTIONS FOR USE

- Plug the appliance to the power supply and slide the function switch to the desired position:

0 – Off.

1 – VIBRATION MASSAGE

- Enclosed attachments allow providing the right massage mode for each part of your body.
- **NOTE:** Although the massager may be applied through clothing, it is most effective when acts directly on the body.

2 – INFRA-RED HEATING

- Infrared radiation penetrates deep into muscles, enhancing blood circulation and reducing tension in muscle areas. Use the infrared heat function can also accelerate the restorative processes in the case of colds, inflammation and minor injuries.

3 – VIBRATION MASSAGE WITH INFRA-RED HEATING

- Combination of infrared heat and vibration massage intensifies effect of both functions.

CLEAN AND CARE

- Switch the appliance off and ensure that the massager has completely cooled.
- Wipe over the housing with a soft damp cloth and then with a dry cloth.
- Do not use hard abrasives or cleaners.
- Attachments may be cleaned in soapy water

STORAGE

- Be sure that the appliance is unplugged and cooled down.
- Complete instructions of section CLEAN AND CARE
- Keep the appliance in a cool, dry place.

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