
ZENNET

MASSAGE MATTRESS ZET-780

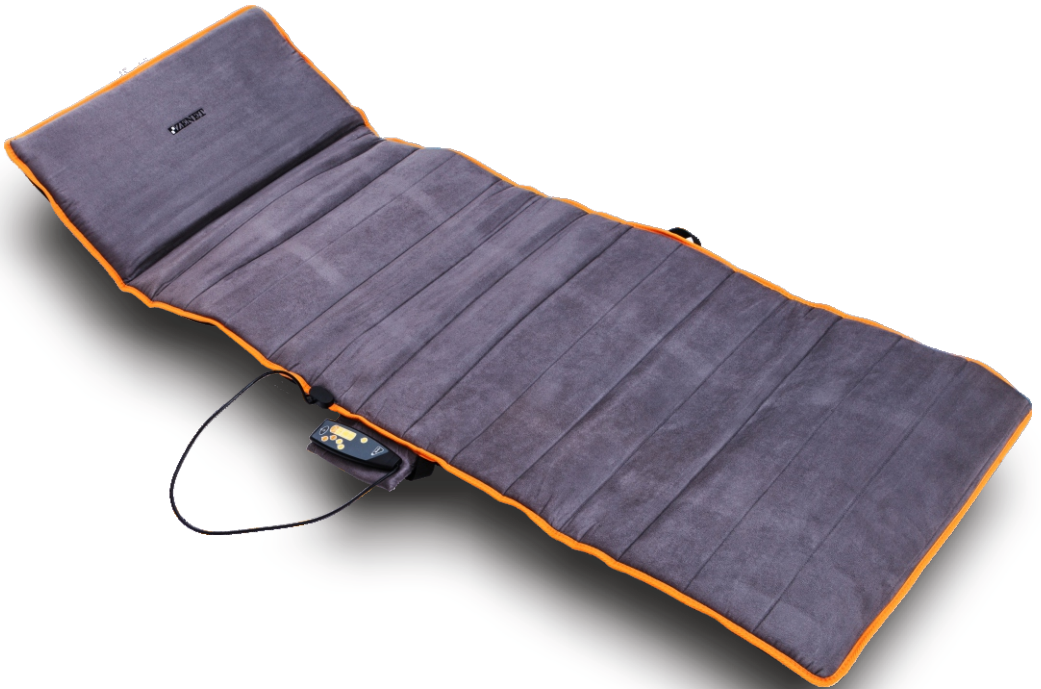MESSAGE MATTRESS ZET-780

INSTRUCTIONS FOR USE

Voltage 220V

Power 14.4W

Dear customer,

Thank you for buying one of
Models of ZENET massagers.

Please read the instruction manual and save it
for future use.

DESCRIPTION:

The massage mattress is equipped with massage vibromotors. The healing effect is greatly enhanced by the heating function. In combination with the warming of the movement of the vibromotors, they relax the muscles and eliminate their cramps. This is an additional way to stimulate blood circulation and the lymphatic system.

CHARACTERISTICS:

1. A pillow in the area of the head of the mattress for comfort.
2. The mattress is made of soft and comfortable fabric.
3. Heat function during the massage.
4. There are five massage programs.
5. Eight massage heads - two for each zone.
6. A relaxing vibration massage for the whole body.

EXPLOITATION:

1. Connect the adapter with the device to a power source.
2. The «POWER» key switches the massager on and off.
3. Head and neck massage - the «NECK» button
4. The «BACK» key activates the back massage function.
5. The «THIGHS» button activates the massage function of the hips and buttocks.
6. The «HEAT» key activates the heating function and the corresponding display lights up.

GUARANTEE LIABILITIES:

1. Guarantee of the device according to the operating rules - 1 year from the date of sale.
2. During the warranty period, the device must be replaced if any defects caused by the manufacturer are found. The sales date and the stamp of the trade organization are required.
3. The guarantee does not apply to devices with:
 - mechanical damage;
 - defects resulting from a violation of the terms of use or incorrect actions by the owner of the device;
 - deficiencies due to natural disasters;
 - defects caused by the ingress of foreign bodies and liquids as a result of self-repair or self-introduction of structural changes;
 - used for commercial or industrial purposes.

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